





## ADVANCED DETECTION METAL DETECTIBLE GASKETS

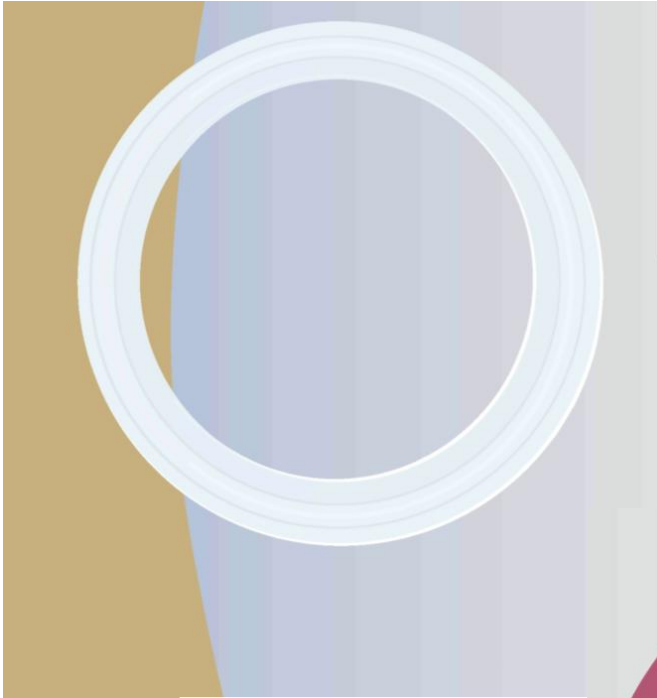

CSKBIO supplies Advanced Detection Gaskets using pedigree of polymers including:

- Platinum Silicone
- PolySteel (Tefon®)
- BUNA
- EPDM
- Viton®

The Polymers are offered in various colors. Please inquire for standard and custom colors.

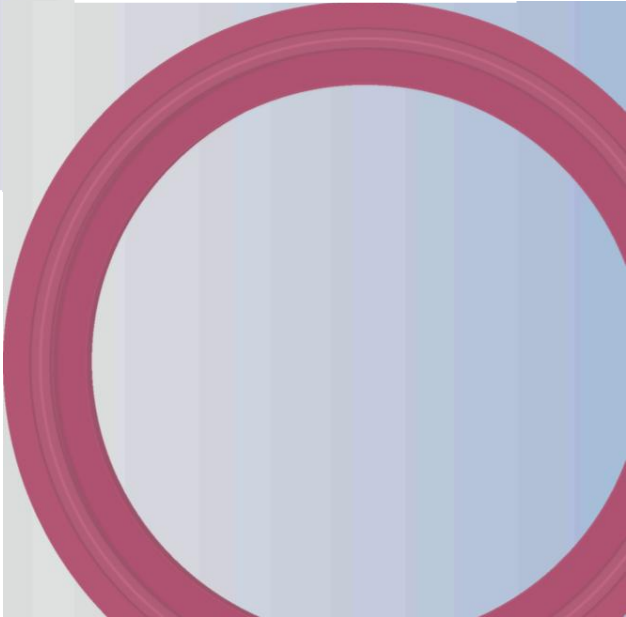



Advanced Detection Gaskets automatically detect polymer decomposition inside your process line.

Our Advanced Detection Program is larger than ever.

We now offer:

AS-568 O-RING

Tri-Clamp

Gaskets

Flanged Gaskets

Bevel Seat John Perry

Q-Lines

I-Line Gaskets

We can also make custom geometries.



CSKBIO's Advanced Detection Gasket Program (ADG) identifies polymer particles that are in your product flow by signaling through standard metal detectors. So users are notified in "ADVANCE" of shipping packaged product to retail stores, diminishing costly and embarrassing product recall. With ADG Processors, we are able to locate gasket particles in-line using metal detectors. X-raying product is no longer necessary to find contaminants.

**CSK-BIO**

Internet: [www.cskbio.com](http://www.cskbio.com)

Email: [sales@cskbio.com](mailto:sales@cskbio.com)



## Considerations when buying Advanced Detection Gaskets:

### Materials:



**Buna**



**EPDM**



**Viton®**



**Platinum Silicone**



**Teflon®/Polysteel**

Teflon® and Viton® are registered trademarks of E.I. du Pont de Nemours and Company

### Color:

**Black Red Blue White Green Mustard Purple**


### Type and Sizes:




**Clamps**  
½" thru 12"



**Flanged Clamps**  
1" thru 12"

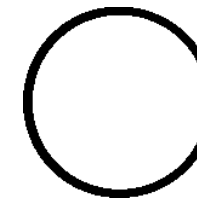


**I-Line**  
1" thru 4"

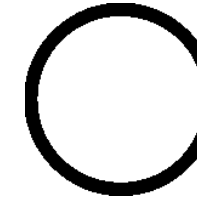

**Bevel Seat**  
1" thru 4"



**John R. Perry**  
1" thru 4"

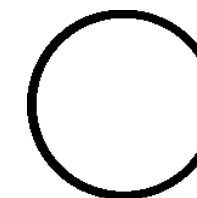

**AS-568 O-Rings**  
EVERY SIZE  
For Pull Pumps, Valves  
and assemblies o all types



**Q-Line**  
1" thru 4"

**CSK-BIO**

Internet: [www.cskbio.com](http://www.cskbio.com)


Email: [sales@cskbio.com](mailto:sales@cskbio.com)



## Advanced Detection Gasket Part Numbers:

\*Advanced Detection Gaskets start with "CSK40"

Gasket Typ	Material Code	Gasket Size	Gasket Codes
<b>CSK40 MP</b> Clamp	<b>G-PS</b> PolySteel	<b>050</b> 1/2"	<b>400</b> 4" <b>BK</b> Black
<b>CSK40 MPF</b> Flanged	<b>U</b> Buna	<b>075</b> 3/4"	<b>500</b> 5" <b>RD</b> Red
<b>CSK40 BS</b> Bevel Seat	<b>E</b> EPDM	<b>100</b> 1"	<b>600</b> 6" <b>BE</b> Blue
<b>CSK40 JP</b> John Perry	<b>SFY</b> Viton AR	<b>150</b> 1.5"	<b>800</b> 8" <b>PE</b> Purple
<b>CSK40 Q</b> Q-Line	<b>VGF</b> Viton GF-600SR	<b>200</b> 2"	<b>1000</b> 10" <b>WE</b> White
<b>CSK40 I</b> I-Line	<b>X</b> Silicone	<b>250</b> 2.5"	<b>1200</b> 12" <b>MD</b> Mustard
		<b>300</b> 3"	

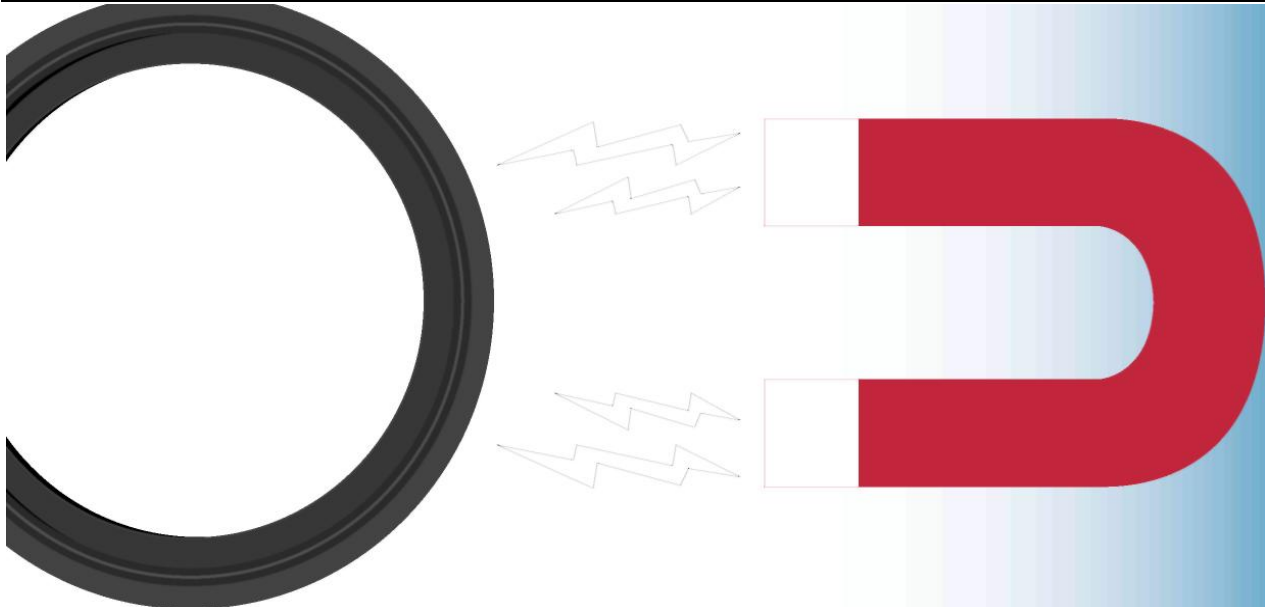

### Examples:

<b>CSK40</b>	MP	U	200	RD	Advanced Detection Clamp Gasket Buna 2"
<b>CSK40</b>	I	E	150	BK	Red Advanced Detection I-Line Gasket EPDM 1.5" Black
<b>CSK40</b>	MPF	SFY	200	BE	Advanced Detection Flanged Gasket Viton R 2" Blue
<b>CSK40</b>	MP	G-PS	50		Advanced Detection PolySteel Gasket